


三環系抗うつ薬

Tricyclic antidepressants

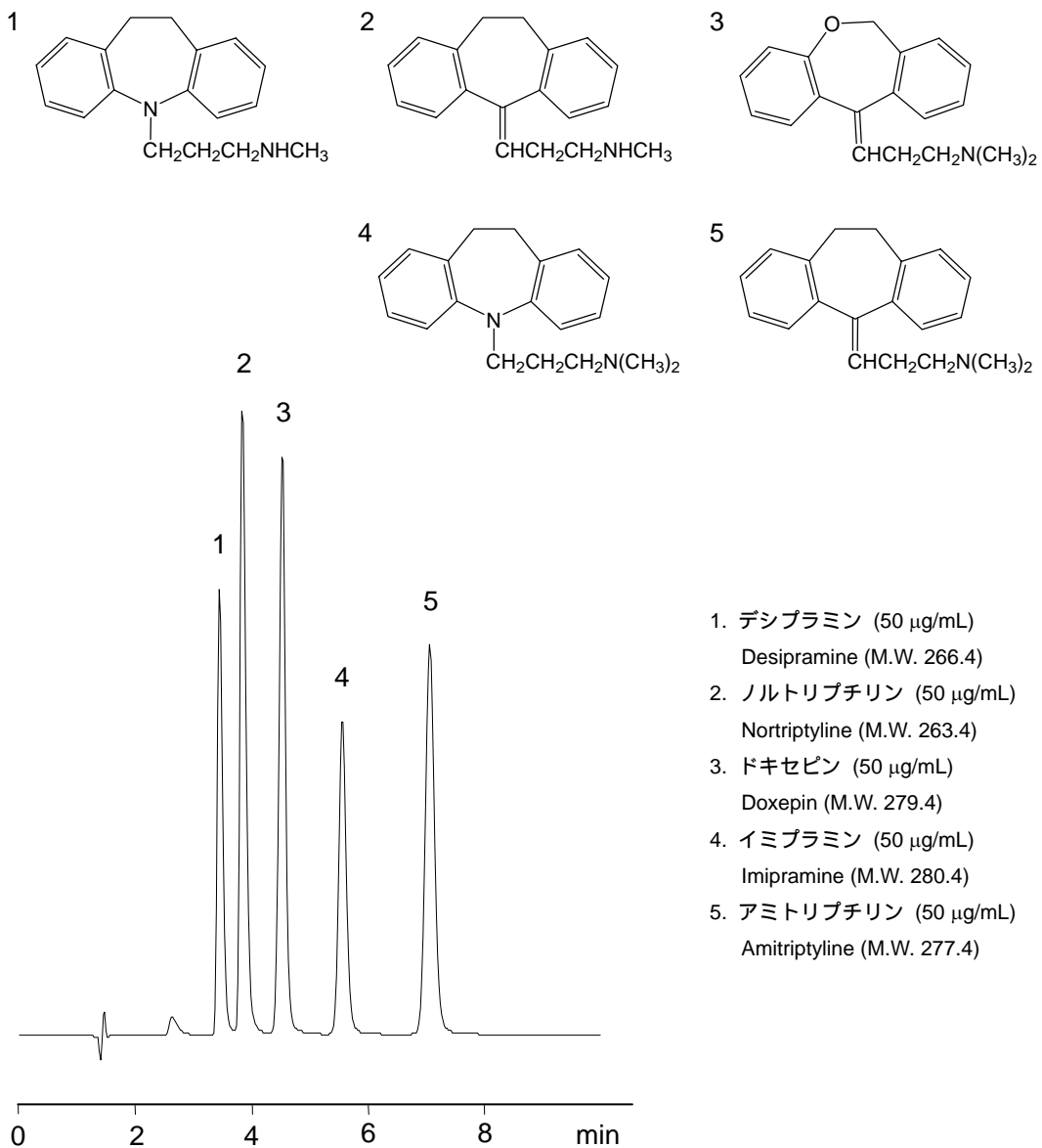

[HPLC Conditions]

Column : CAPCELL PAK C₁₈ MGII S5 ; 4.6 mm i.d. x 150 mm
 Mobile phase : 10 mmol/L CH₃COONH₄ / CH₃OH = 30 / 70
 Flow rate : 1 mL/min
 Temperature : 40
 Detection : UV 210 nm
 Inj. vol. : 1 μL
 Sample dissolved in : CH₃OH
 1 μg/mL = 1 ppm