
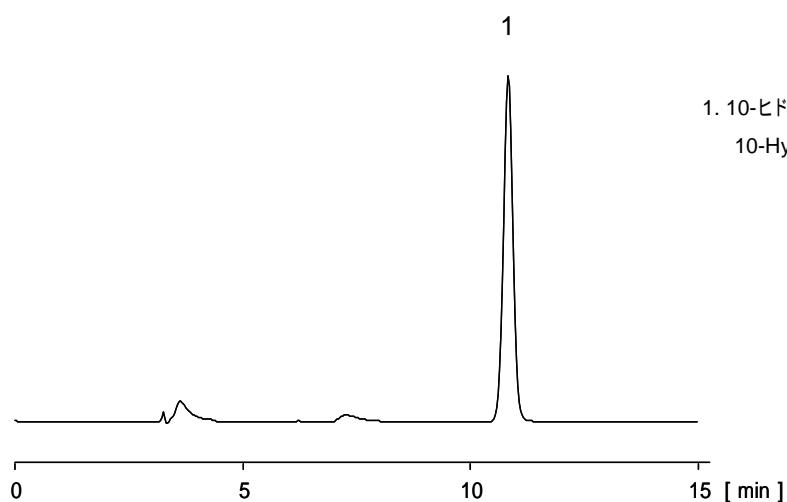


10-ヒドロキシ-2-デセン酸

10-Hydroxy-2-decenoic acid

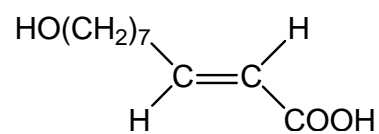

1. 10-ヒドロキシ-2-デセン酸 (100 µg/mL)  
10-Hydroxy-2-decenoic acid (M.W. 186.3)

**[HPLC Conditions]**

Column : CAPCELL PAK C<sub>18</sub> UG120 S5 ; 2.0 mm i.d. x 250 mm  
Mobile phase : 1.5 vol% H<sub>3</sub>PO<sub>4</sub> / CH<sub>3</sub>OH = 50 / 50  
Flow rate : 200 µL/min  
Temperature : 40 °C  
Detection : UV 210 nm  
Inj. vol. : 2 µL  
Sample dissolved in : Mobile phase  
1 µg/mL = 1 ppm