


Sodium alkyl sulfonate

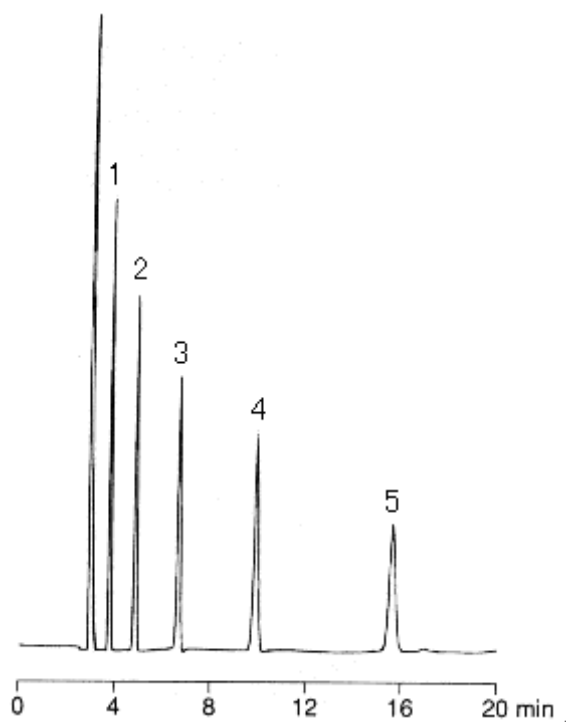

1	Sodium decyl sulfonate (C10) 1.0 %	2	Sodium lauryl sulfonate (C12) 1.0 %	3	Sodium myristyl sulfonate (C14) 1.0 %
4	Sodium cetyl sulfonate (C16) 1.0 %	5	Sodium stearyl sulfonate (C18) 1.0 %		

HPLC Conditions

Column	CAPCELL PAK C18 AG120 S5 : 4.6 mm i.d.×250 mm
Mobile phase	1mol/L NaClO ₄ , CH ₃ OH / H ₂ O = 90 / 10
Flow rate	1.0 mL/min
Temperature	50 °C
Detect	RI
Inj. vol.	10 μL