


Data provided by : Dr.Oka, Aichi Prefectural Institute of Public Health

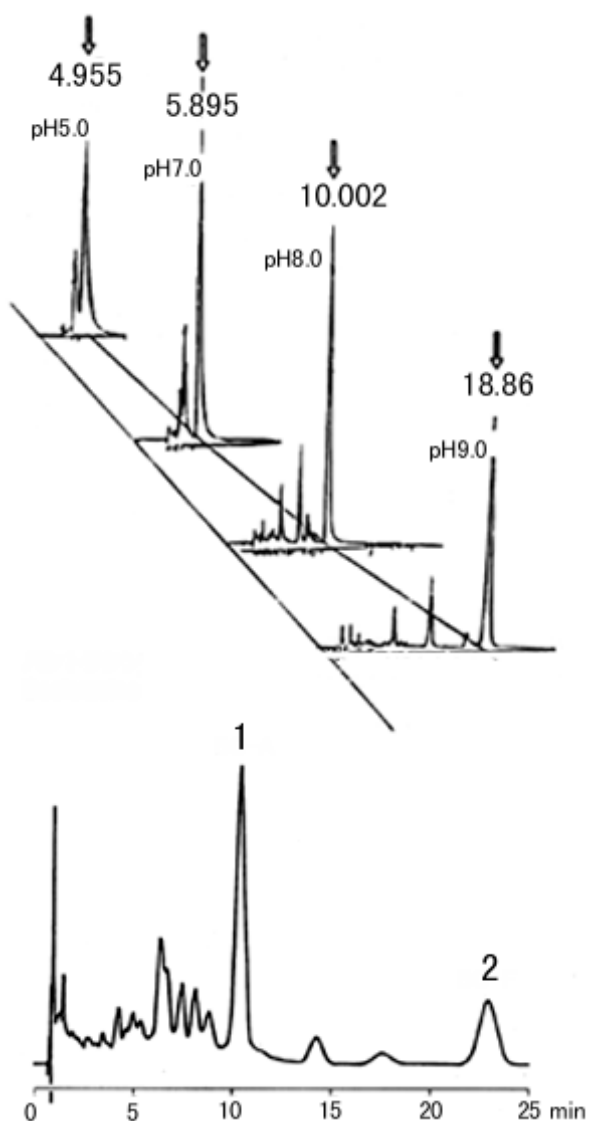

1 Bacitracin-A

2 Bacitracin-F

HPLC Conditions 1

Column CAPCELL PAK C₁₈ AG120 S5 : 4.6 mm i.d.×250 mm

Mobile phase 50mmol/L NaH₂PO₄ / CH₃CN

Detect UV 210 nm

HPLC Conditions 2

Column CAPCELL PAK C₁₈ AG120 S5 : 4.6 mm i.d.×250 mm

Mobile phase CH₃OH / 40mmol/L Na₂HPO₄ = 60 / 40, pH9-10

Detect

UV 234 nm