


Date provided by: Dr.Samejima, Pharmacy Hoshi University

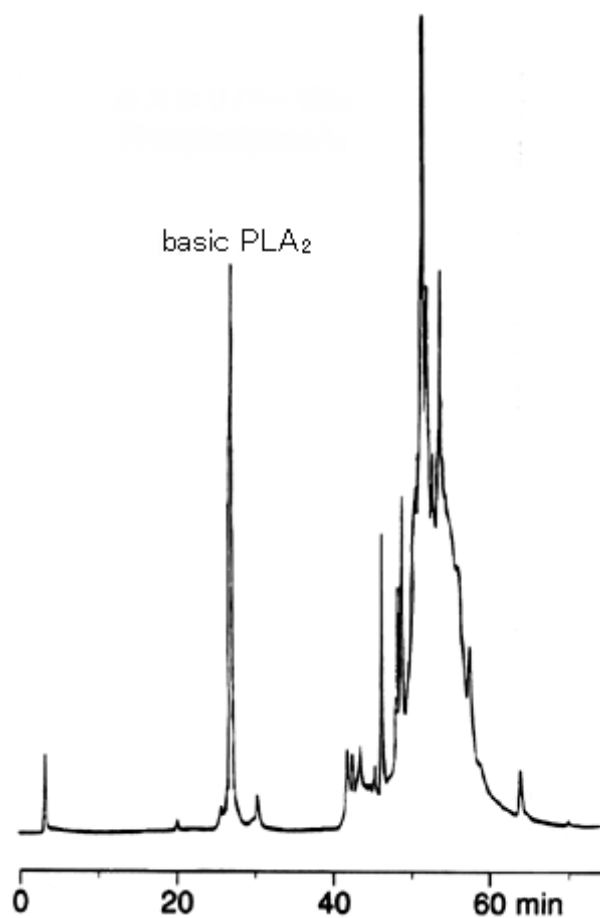

## 1 PhospholipaseA2

### HPLC Conditions

Column	CAPCELL PAK C18 SG300 S5 : 4.6 mm i.d.×250 mm
Mobile phase	(A) 0.1vol% TFA, H <sub>2</sub> O (B) 0.1vol% TFA, CH <sub>3</sub> CN B 25% - 35% (30min) - 90% (40min) Gradient
Flow rate	1.0 mL/min
Detect	UV 280 nm