
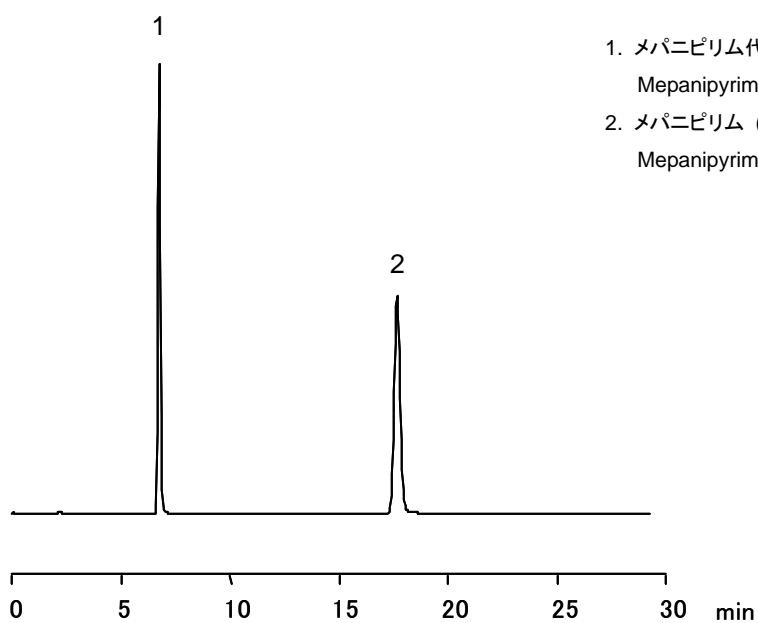
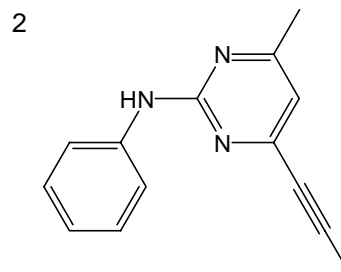


## メパニピリム

## Mepanipyrim

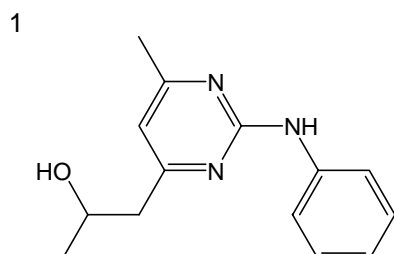

1. メパニピリム代謝物 (50 ppm)  
Mepanipyrim metabolite (M.W. 243.3)
2. メパニピリム (50 ppm)  
Mepanipyrim (M.W. 223.3)

## 【HPLC Conditions】

Column : CAPCELL PAK C<sub>18</sub> MGII S5 ; 2.0 mm i.d. x 250 mm  
Mobile phase : H<sub>2</sub>O / CH<sub>3</sub>CN = 1 / 1  
Flow rate : 200 μL/min  
Temperature : 40 °C  
Detection : UV 270 nm  
Inj. vol. : 2 μL  
Sample dissolved in : CH<sub>3</sub>CN (50 ppm)  
※ 1 μg/mL = 1 ppm