


非ステロイド系抗炎症薬

Nonsteroidal anti-inflammatory drugs

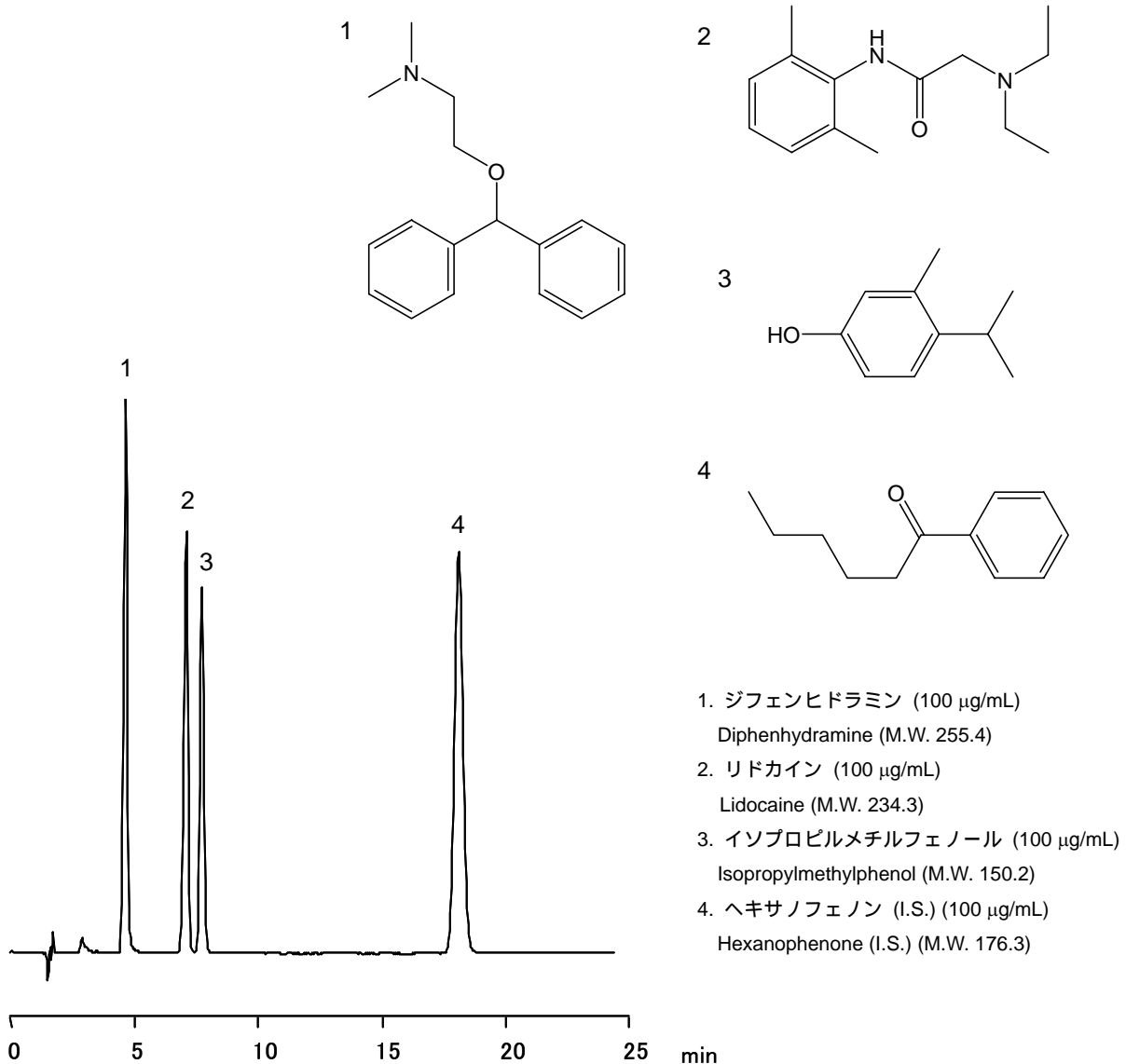

【HPLC Conditions】

Column : CAPCELL PAK C₁₈ MGII S5 ; 2.0 mm i.d. x 150 mm
 Mobile phase : 10 mmol/L CH₃COONH₄ / CH₃OH = 40 / 60
 Flow rate : 200 μL/min
 Temperature : 40 °C
 Detection : UV 230 nm
 Inj. vol. : 2 μL
 Sample dissolved in : 50 vol% CH₃CN

※ 1 μg/mL = 1 ppm