


水溶性ビタミン (イオンペア法) Hydrophilic vitamins (ion-pair method)

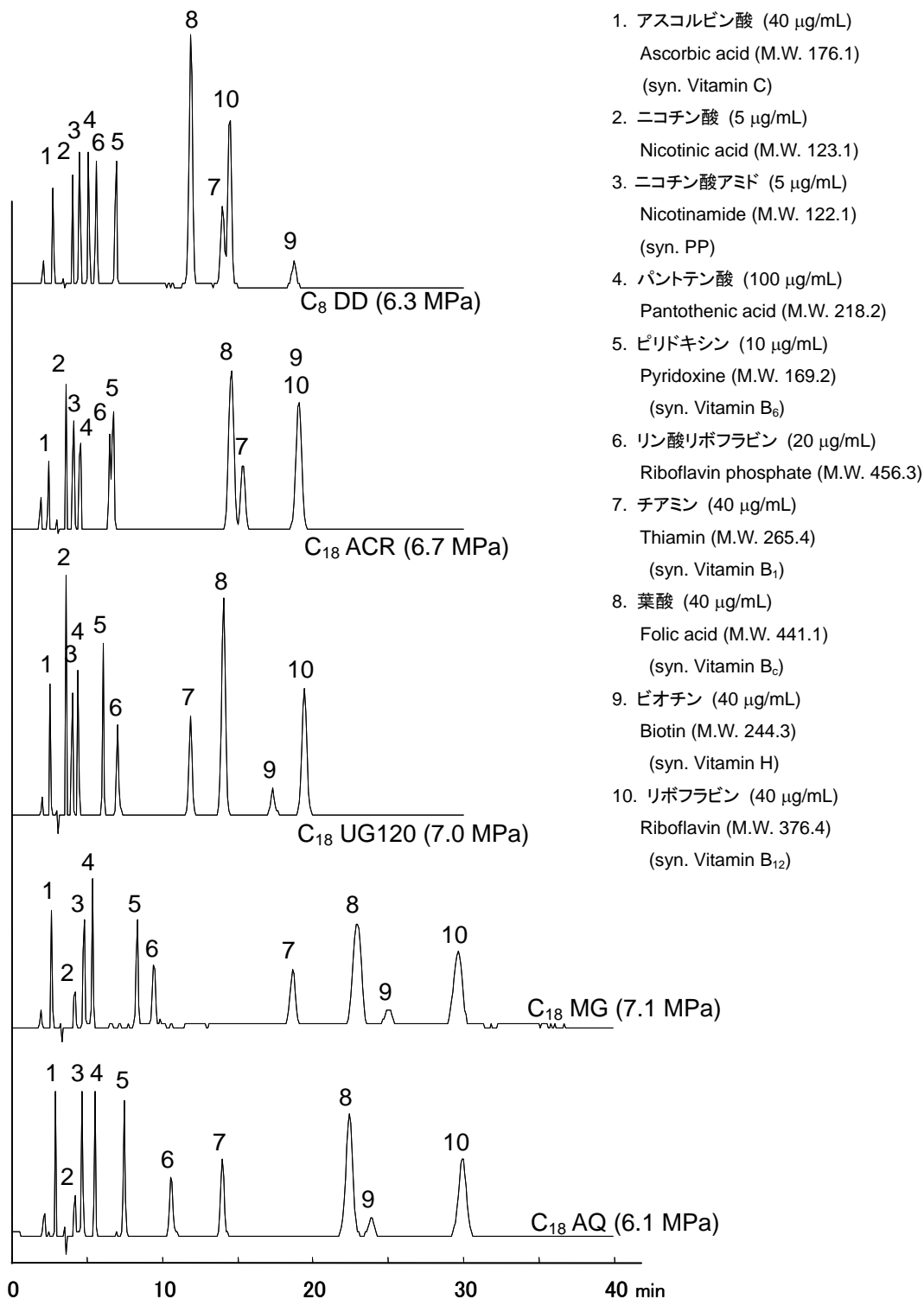

【HPLC Conditions】

Column size : 4.6 mm i.d. x 150 mm
Mobile phase : Sodium hexane sulfonate was added to a solution (20 mmol/L H₃PO₄ / CH₃CN = 91 / 9) at 5 mmol/L.
*The amount of acetonitrile significantly affects retention times in this application.
Flow rate : 1 mL/min
Temperature : 40 °C
Detection : UV 210 nm
Inj. vol. : 10 µL
Sample dissolved in : Refer to the page, "Preparation of hydrophilic vitamin solutions".
(Final concentrations should be changed as needed.)
※ 1 µg/mL = 1 ppm