
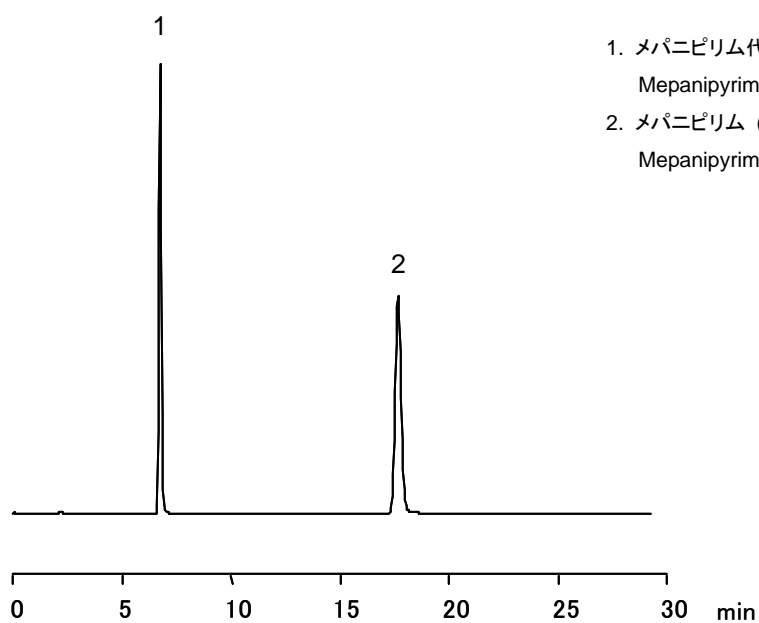
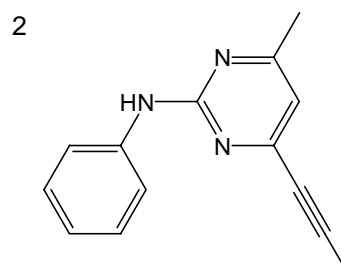


メパニピリム

Mepanipyrim

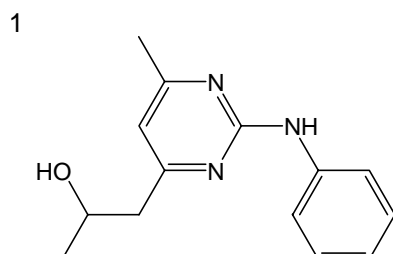

1. メパニピリム代謝物 (50 ppm)
Mepanipyrim metabolite (M.W. 243.3)
2. メパニピリム (50 ppm)
Mepanipyrim (M.W. 223.3)

【HPLC Conditions】

Column : CAPCELL PAK C₁₈ MGII S5 ; 2.0 mm i.d. x 250 mm
Mobile phase : H₂O / CH₃CN = 1 / 1
Flow rate : 200 μL/min
Temperature : 40 °C
Detection : UV 270 nm
Inj. vol. : 2 μL
Sample dissolved in : CH₃CN (50 ppm)
※ 1 μg/mL = 1 ppm